



D-U2000
Flexible Veterinary Dental Unit

RWD



**EASY
FLEXIBLE**

RWD Veterinary Dentistry , Easy & Safe

D-U2000

Flexible Veterinary Dental Unit

D-U2000 Flexible veterinary dental unit was designed to help veterinarians perform their daily work easier and more efficient. Veterinarians can complete grinding, polishing, cutting and cleaning very quickly without fatigue.

With RWD dental unit, you will...

- ① Get high return on equipment investment.
- ② Shorten the treatment, therefore reduce the risk of anesthesia.
- ③ Reduce the damage of alveolar bone and oral mucosa caused by violent tooth extraction.
- ④ Reduce the risk of retained roots.
- ⑤ Make more flexible dental treatment solutions

28000Hz scaler with **automatic frequency tracking** to ensure power, and **white LED light** for clear vision without fatigue.

Low speed motor with coolant delivery tube for **effective cooling**. It can also be compatible with polishing head and burs.


Whole machine **weighs only 32kg**. Super smooth casters will ensure **flexible movement**.

Magnetic platform to absorb burs or other fine tools, **avoid accidental sliding or loss**.


Powerful 360° swivel high speed handpiece with **self-generating LED illumination** helps veterinarians cut animal teeth easily

Automatically drain pollution or water from the air tank, greatly **reduces the trouble of daily maintenance for veterinary**, reduces the risk of rust, and prolong the lifespan



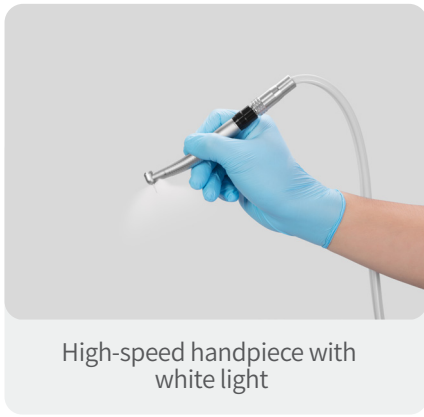

High-speed handpiece with white light

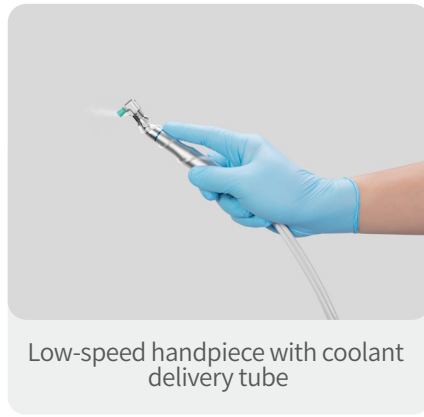

Low-speed handpiece with coolant delivery tube




Operation console



3 way air/ water syringe



Magnetic platform

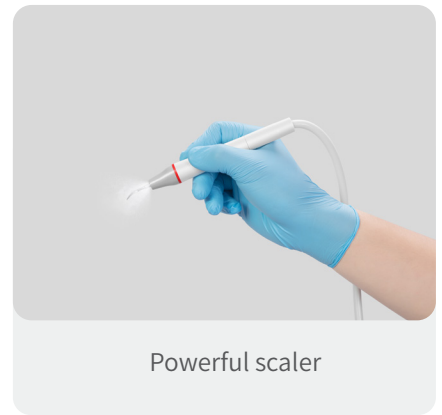

Powerful scaler

Technical parameters

Item	Specification	
Dimension	Length	45cm
	Width	40cm
	Height	80-105cm
	Weight	32KG
Oil-Free Air Compressor	Power	110/220V, 50/60Hz
	Out Put Pressure	0.45-0.5Mpa
High Speed Handpiece	330000-390000rpm	
Low Speed Handpiece	5000-20000rpm	
Auto Drain Function	Single time	0.5~10s
	Interval time	0.5~45min
Sound Level	≤65dB	

Main configuration

Description	D-U2000 (Standard Edition)	D-U2000L (Lite)
High Speed Handpiece	1pcs	—
Low Speed Handpiece	1pcs	—
3 Way Air/Water Syringe	1pcs	1pcs
Scaler Handpiece	1pcs	1pcs
Burs	1 Set (Tapered fissure burs*5+Spherical burs*5)	—
Disposable Polishing Head	1 Package (5pcs)	—
Scaling Tips	1 Set (5pcs)	1 Set (5pcs)

Note: The D-U2000L (Lite) reserves the interface for high and low speed handpiece

One Stop Solution!



Anesthesia System



Veterinary Dental Unit



Veterinary Dental X-Ray

RWD Life Science Co.,Ltd

Add: 19-20/F, Building 7A, Shenzhen International Innovation Valley, Dashi 1st Road, Nanshan District, Shenzhen, Guangdong, P.R. China E-mail: market@rwdstco.com Web: www.rwdstco.com

RWD Life Science Inc.

Add: 850 New Burton Road, Suite 201, Dover, DE 19904 (New Office)

Tel: (858)900-5879 E-mail: service@rwdls.com Facebook:@RWDveterinary

Linkedin: RWD Life Science Twitter:@RWDScience

RWD Flexible Veterinary Dental Unit-V1.0-201102-EN